


Product Datasheet

Goat anti-human CRP antibody

#SB-Pab-0304

Immunogen	Purified human C reactive protein
Quantity	0.2mg in 1mg/ml
Purification	Ammonium sulphate precipitation of serum followed by antigen affinity purification

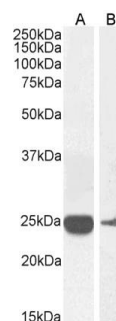

Buffer	Phosphate buffer saline, PH 7.2, with 0.02% sodium azide as preservative	<i>SB-Pab-0304</i> anti-CRP (0.3µg/ml) staining of Human Liver (A) and Human Serum (B) lysates (35µg protein in RIPA buffer). Primary incubation was 1 hour. Detected by chemiluminescence
Storage and handling	Store at 4C for short durations. Store at -20 for longer periods (> 1 month). Minimize freeze thaw cycles	
Applications tested	Protein ELISA: antibody detection at 1:128000 dilution Western blotting : ca 25kDa band in Human Liver and Serum lysates	
Application notes	Please optimize dilutions for your particular assay conditions. Starting test dilution: 0.5µg/ml for ELISA and WB	