


Product Datasheet

Goat anti-human CRP antibody-biotin

#SB-Pab-0304B

Immunogen	Purified human C reactive protein
Quantity	0.2mg in 0.5mg/ml
Purification	Ammonium sulphate precipitation of serum followed by antigen affinity purification

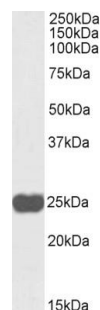

Buffer	Phosphate buffer saline, PH 7.2, with 0.02% sodium azide as preservative
Storage and handling	Store at 4C for short durations. Store at -20 for longer periods (> 1 month). Minimize freeze thaw cycles

SB-Pab-0304B (0.3µg/ml) staining of Human Liver lysate (35µg protein in RIPA buffer). Primary incubation was 1 hour. Detected by chemiluminescence

Applications tested	Protein ELISA: antibody detection at 1:128000 dilution Western blotting : ca 25kDa band in Human Liver lysate
----------------------------	--

Application notes	Please optimize dilutions for your particular assay conditions. Starting test dilution: 0.5µg/ml for ELISA and WB
--------------------------	---